
THOR X CORE

SENIOR TECH

*Static Pedal Exerciser / Pedaleador estático / Pedale
static / Statischer Pedaltrainer / Pédalier statique*

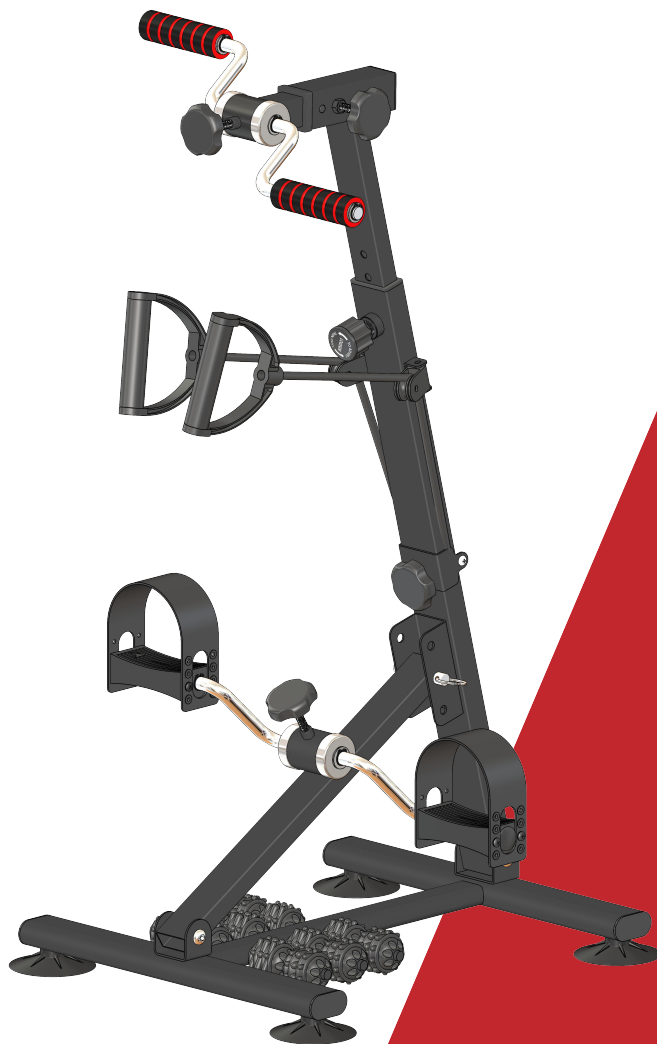INSTRUCTIONS /
INSTRUCCIONES /
ISTRUZIONI /
ANWEISUNGEN /
INSTRUCTIONS

Product description

The THORXCORE SENIOR TECH Static Pedal Exerciser is designed for gentle home exercise of the arms and legs. It promotes circulation, mobility, and muscle activity. This product is not a medical or rehabilitation device and is intended for home use by adults only.

Safety warnings

- Use only on a flat, stable surface.
- Keep children and pets away. This product is not a toy
- Ensure all screws, knobs, and straps are securely tightened before each use.
- Keep fingers, clothing, and hair away from moving parts.
- Adjust resistance gradually. Stop exercising immediately if you feel pain, dizziness, or discomfort.
- Do not use the product if any part is damaged, loose, or missing.
- This product is not a medical or rehabilitation device. Consult your doctor before use if you have any health condition.
- Place the exerciser on a non-slip surface and use indoors only.
- Do not modify or disassemble the product.
- Keep this manual for future reference.

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Parts and accesories included:

Pictures	Part Number	Components	Quantities	Pictures	Part Number	Components	Quantities
	A	Base frame	1		H	Seven-star Plum Blossom Knob	1
	B	Connecting Bracket	1		I	Latch	1
	C	Top bar	1		J	Screw	1
	D	Hand bracket	1		K	Lock nut	1
	E	Non-slip Foot Mats	4		L	Open-End Wrench	1
	F	Pull-out Knob	1		M	L-Shape Wrench	1
	G	Plum Blossom Knob	3				

1

EN: Place all parts on a clean, flat surface /

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2

EN: Attach the 4 Foot mats (E) to the Base frame (A) and tighten with the wrench (L) if needed /

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3

EN: Insert one plum knob (G) into the pedal axis. This knob can be adjusted as required in order to add or decrease the resistance of pedals rotation for higher/lower exercise intensity /

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4

EN: Insert the connecting bracket (B) into the base frame and secure with the screw (J) and lock nut (K). /

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5


EN: Insert one plum knob (G) into the corresponding hole of the connecting bracket to keep it extended vertically /

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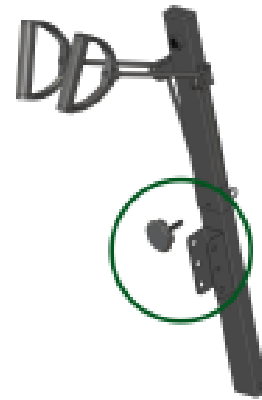6

EN: Attach the pedals bar to the connecting bracket using the Latch (I) provided

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7

EN: Insert the top bar (C) into the connecting bracket and fix with the Pull knob (F). Ensure that the hole of the top bar is headed towards the pedals.

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9


EN: Insert one plum knob (G) into the hand pedal axis. This knob can be adjusted as required in order to add or decrease the resistance of pedals rotation for higher/lower exercise intensity

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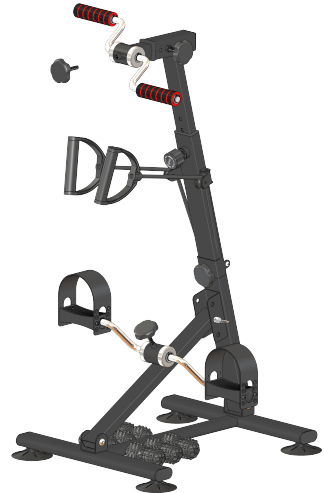10

EN: The hand bracket can also be positioned vertically, allowing the hand pedals to rotate on a different axis for an alternative upper-body exercise. /

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ENGLISH

Use

- Place the device on a stable surface, whether standing or seated.
- For legs: pedal slowly and increase resistance as desired.
- For arms and shoulders: hold the upper handles and move them gently back and forth.
- Use the elastic bands for upper-body exercises. Do not overstretch them.
- Recommended duration: start with 10–15 minutes per session and increase gradually.

Maintenance

- Clean with a soft, slightly damp cloth. Do not use solvents or immerse in water.
- Check screws, straps, and rubber caps periodically.
- Store in a dry place, away from heat and humidity.
- Replace any damaged or worn parts immediately.

Disposal and Legal Notice

Dispose of the packaging in the cardboard recycling container.

Do not discard with household waste.

This product is intended for home exercise only. Keep this manual for future reference.

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